


# Switch Adapted Bubble Blower

## USER QUICK GUIDE

### Introduction

The Switch Adapted Gazillion Tornado Bubble Machine is an adapted toy that has a 3.5 mm jack. To activate the device, an assistive switch with a 3.5 mm plug can be plugged into the jack on the bubble blower. The Switch Adapted Bubble Blower may be beneficial to someone who may have difficulty pressing the original small activation button due to its size, location, or required activation force.

### Features

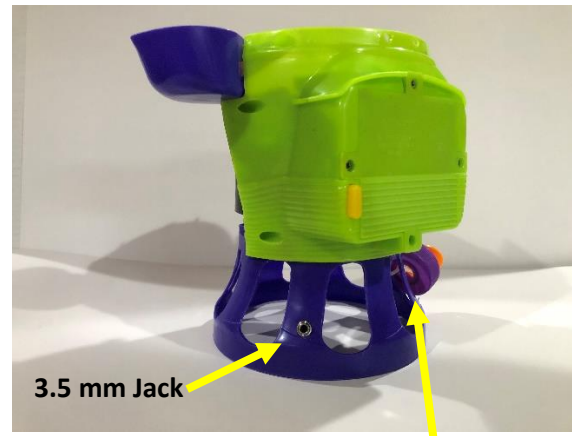

**Bubble cover with hole in middle**

### Usage

**Step 1:** Remove cap from the bottle of bubbles.

**Step 2:** Put the alternative purple and orange cover with the hole in the middle on the bottle of bubbles.

**Step 3:** Flip the bottle of bubbles 180 degrees and insert the bottle of bubbles into the purple part of the bubble machine as shown below.



# Switch Adapted Bubble Blower

## USER QUICK GUIDE

**Step 4:** Plug an assistive switch into the 3.5 mm jack on the back side of the bubble machine.

**Step 5:** Activate the assistive switch or original switch to generate bubbles. It may take a few seconds before bubbles are generated. The assistive switch will need to be held down to continue to generate bubbles.

### Compatibility

The bubble machine should work with any assistive switch that has a 3.5 mm plug.

### Specifications

Length (mm)	Width (mm)	Height (mm)	Mass (g)
165	165	150	725

### Cleaning

The device can be wiped down with a damp cloth. Any remaining bubbles in the device can be poured back into a bottle for reuse.